



**Don't be Puzzled by Glaucoma –  
The Comprehensive Solution for Your Practice**

**HEIDELBERG  
ENGINEERING**



## All Pieces of the Puzzle

### The Complete Picture of Glaucoma Diagnosis

Glaucoma is one of the leading causes of blindness. Early glaucoma detection and reliable progression analysis may prevent loss of vision. A serious management of glaucoma is the basis of an individualized therapy and excellent patient care.

The detection and monitoring of glaucoma require a targeted combination of clinical observations and reliable examination results. Clinical expertise is always the key to an accurate glaucoma diagnosis.

Structural and functional examinations with Heidelberg Engineering products provide invaluable information to aide in the diagnosis and treatment of glaucoma.

Advanced glaucoma diagnosis demands advanced technology. A complete picture of glaucoma develops only if the diagnostic information fits together like the pieces of a puzzle.

We provide serious diagnostic power for glaucoma!  
Comprehensive glaucoma solutions from Heidelberg Engineering.









## Flexible Integration – Efficient Networking

### Intelligent Use of Patient Data and Optimized Patient Workflow

Efficient integration of diagnostic equipment in the clinical setting can significantly improve workflow and management of an increasing number of patients while delivering a high quality of care.


The HEYEX software platform provides flexible solutions from a solo practice to a large clinic environment. The integrated patient database facilitates storage and management of image files from all existing Heidelberg Engineering devices, such as HRT, HEP and SPECTRALIS.

Network viewing software allows easy and reliable access to patient image files on any networked computer supporting seamless consultation and diagnosis.

-  One software platform for all Heidelberg Engineering devices
-  Network integration and viewing stations
-  Data security
-  DICOM compatibility
-  Connectivity to many electronic health records systems
-  Import of third party image data

# Validated Diagnostics

The excellent diagnostic performance of the Heidelberg Retina Tomograph (HRT) is well established. The HRT has been used worldwide over 20 years for the detection and management of glaucoma. Advanced algorithms like the Moorfields Regression Analysis (MRA) and the Glaucoma Probability Score (GPS) allow an objective assessment of structural glaucomatous damage. Objective structural parameters complement the 3D topographic profile of the optic nerve head and facilitate a comprehensive risk assessment of suspect glaucoma patients.

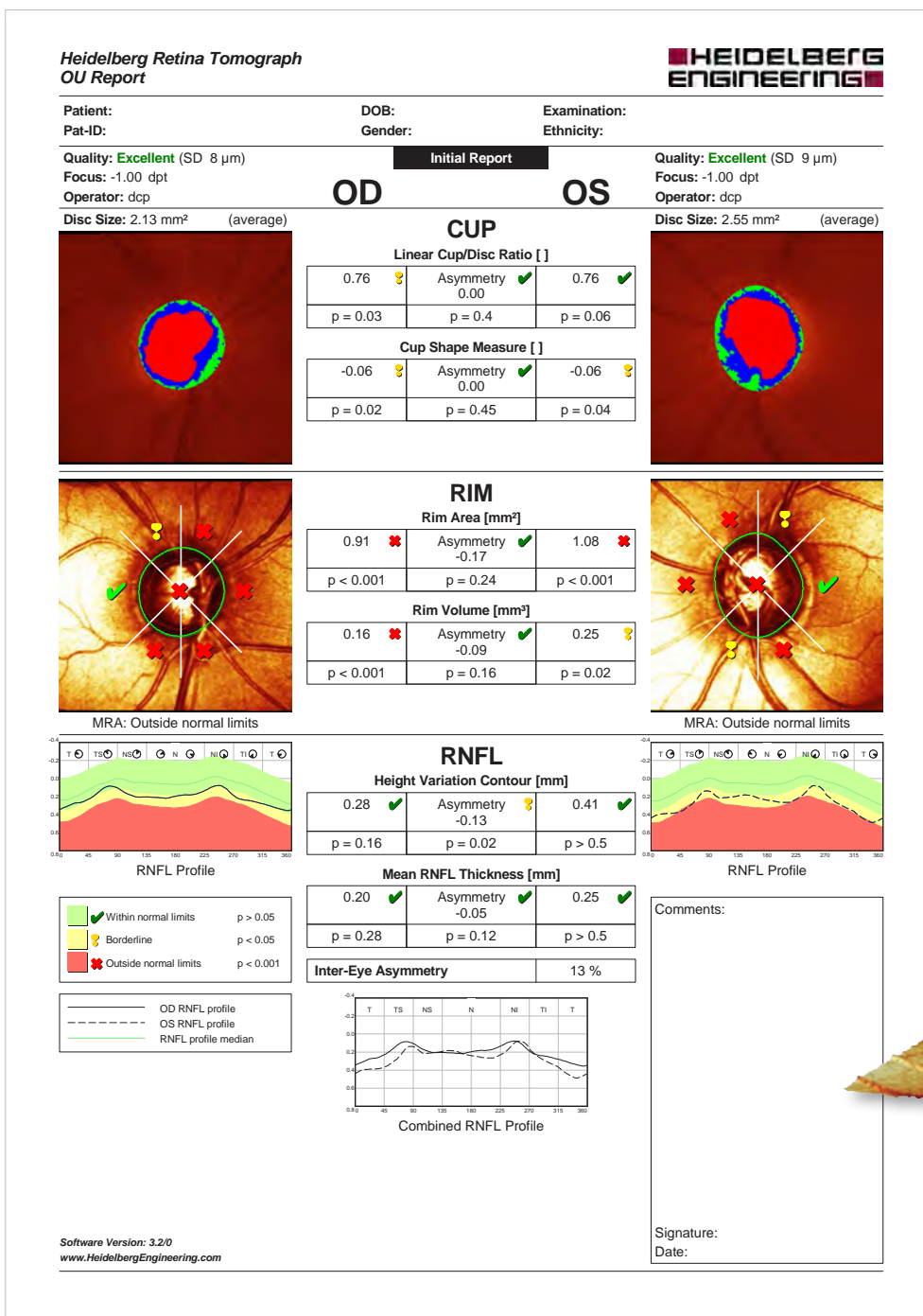


All relevant anatomic structures at a glance – optic nerve head (ONH) excavation, neuroretinal rim and peripapillary retinal nerve fiber layer

Weinreb 2010: Baseline GPS, MRA, and stereoparameters alone or when combined with baseline clinical and demographic factors can be used to predict the development of POAG end points in OHTS (Ocular Hypertensive Treatment Study) participants and are as effective as stereo photographs for estimating the risk of developing POAG in ocular hypertensive subjects<sup>1</sup>

Zangwill 2005: HRT MRA is the top predictive factor for glaucoma and can identify those patients at high risk for developing the disease<sup>2</sup>

Strouthidis 2010: In clinical practice, patients with both an abnormal MRA and an abnormal GPS classification at presentation should be followed closely as they have an increased risk of future visual field or HRT change<sup>3</sup>

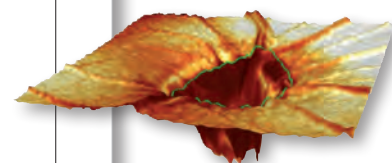

<sup>1</sup> Weinreb RN et al; Ophthalmol 2010; 117:1674–1683.

<sup>2</sup> Zangwill LM. et al.; Arch Ophthalmol 2005;123:1188-1197.

<sup>3</sup> Strouthidis NG et al.; Glaucoma 2010; 19(5):304-309.

## Early Detection and Reliable Progression Analysis

The Topographic Change Analysis (TCA) of the HRT is the best known and most evidence-based method for glaucoma management<sup>4</sup>. The Gold Standard for objective comprehensive structural analysis of the ONH is the key to an individualized approach to glaucoma diagnosis and treatment.


**Heidelberg Retina Tomograph  
TCA Overview** **HEIDELBERG  
ENGINEERING**

---


**Patient:** Sex:      DOB:      Pat-ID:      Ethnicity: **OS**

**Examination:** Baseline:      Last Follow-Up:      Elapsed:


---

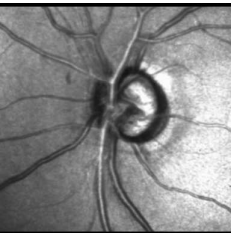


Baseline: 30.Mrz.2005

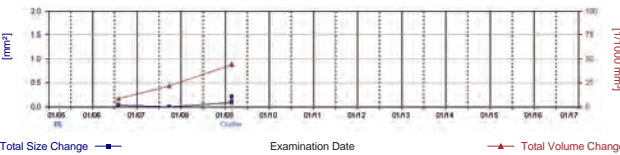


Examination Date

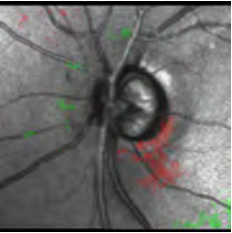

Follow-Up: #2, 01.Aug.2006

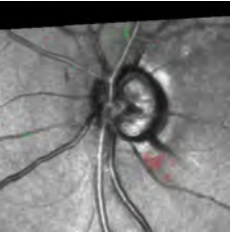

Follow-Up: #3, 26.Sep.2007

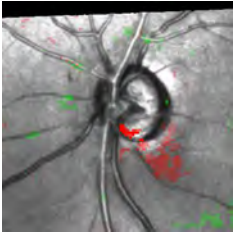



Follow-Up: #4, 23.Feb.2009


- ONH and retinal nerve fiber layer change at a glance
- Early detection of subtle glaucomatous changes i.e. disease progression
- Trend analysis – Volume and area changes indicated
- Contour-line and reference plane independent
- Follow-up examinations and single scans automatically aligned

<sup>4</sup> Chauhan BC. et al.; Invest Ophthalmol Vis Sci 2000; 41:775-782.

# One Perimeter – Two Methods


The visual field examination remains an essential part of glaucoma diagnosis and management, despite rapid developments in imaging technologies<sup>5</sup>. Heidelberg Engineering complements the Gold Standard of structural glaucoma diagnosis, the Heidelberg Retina Tomograph (HRT) with an innovative and reliable partner, the Heidelberg Edge Perimeter (HEP). The HEP offers two perimetric methods in one device.

-  Flicker-Defined-Form (FDF) Perimetry for early glaucoma detection in patients at risk and ocular hypertensive patients
-  The Gold Standard, White-on-White Perimetry (SAP) for continuous care of patients with manifest glaucoma
-  Unique Legacy Scale for a direct comparison with HFA examination results
-  The Functional Change Analysis FCA allows monitoring of significant functional changes over time

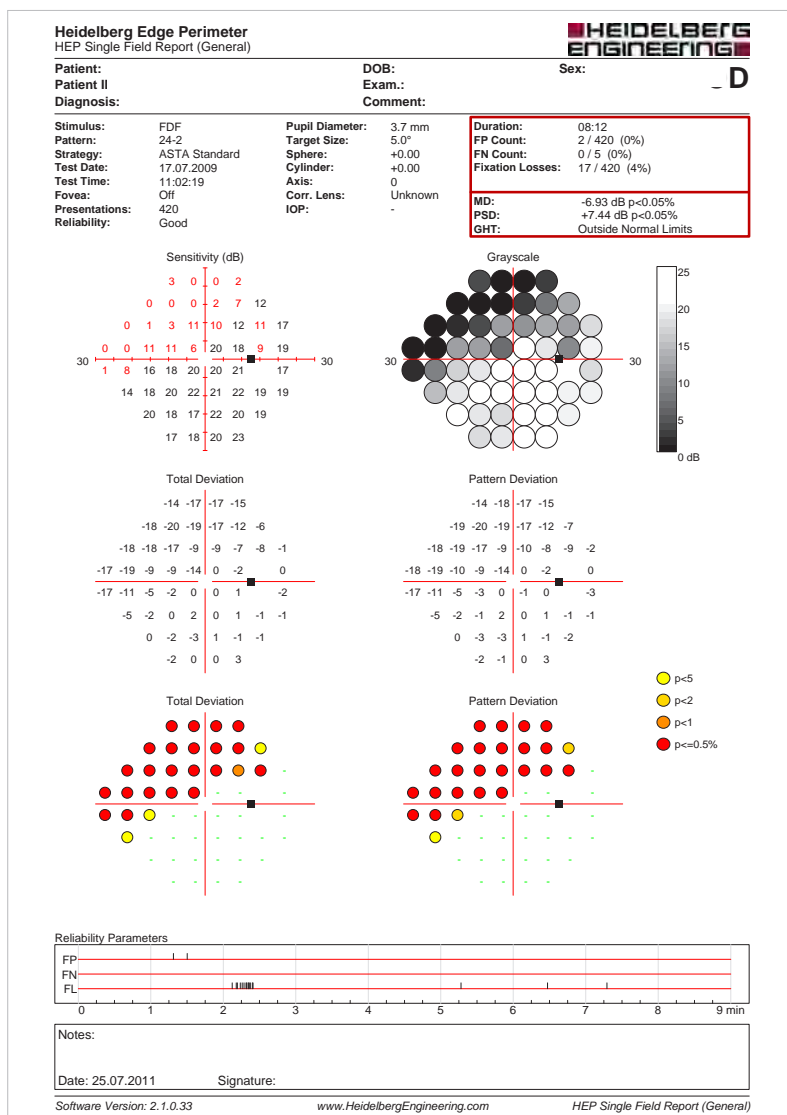

Examination Report – HEP Single Exam Report

Examination Report – Functional Change Analysis FCA

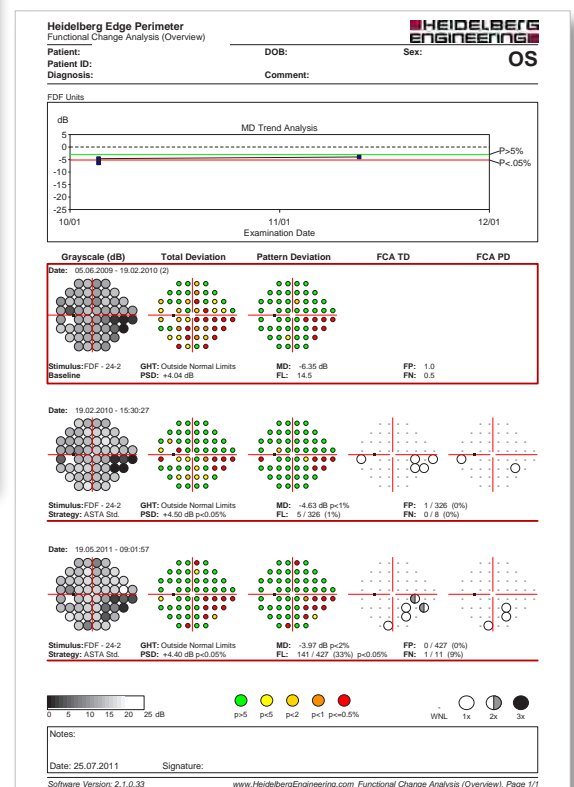

<sup>5</sup> Terminology and Guidelines for Glaucoma. 3rd Edition (2008). European Glaucoma Society. Image courtesy Dr. Fritz Dannheim, Seevetal-Hittfeld, Germany



# Structure and Function

Clinically, glaucomatous damage is typically diagnosed by combining structural and functional assessments. In particular, anatomical structural damage detected by Confocal Scanning Laser Ophthalmoscopy (HRT) or Optical Coherence Tomography (OCT) is compared to changes in the visual field (VF) sensitivity as measured by FDF Perimetry or SAP. Heidelberg Engineering combines datasets of HRT, SPECTRALIS and HEP in the common software platform HEYEX, and conjointly displays structure-function results on a single report.

- All relevant information on one printout
- First real correlation of structural and functional examination results
- Reliable clinical diagnosis
- Improved patient and time management

**HRT + HEP  
OU Report**

**HEIDELBERG  
ENGINEERING**

Patient: \_\_\_\_\_ Sex: \_\_\_\_\_ DOB: \_\_\_\_\_ Pat-ID: \_\_\_\_\_ Ethnicity: \_\_\_\_\_

Examination: HRT Date: \_\_\_\_\_  
HEP(OD) Date: \_\_\_\_\_  
HEP(OS) Date: \_\_\_\_\_

---

**OD**

Disc Size: 2.71 mm<sup>2</sup> (large)

MRA: Borderline

**OS**

Disc Size: 2.30 mm<sup>2</sup> (average)

MRA: Outside normal limits

Test	Value	Reference
Linear Cup/Disc Ratio	0.67 p = 0.24	0.57 p = 0.42
Cup Shape Measure	-0.13 p > 0.5	-0.23 p > 0.5
Rim Area [mm <sup>2</sup> ]	1.48 p = 0.04	1.56 p > 0.5
Rim Volume [mm <sup>3</sup> ]	0.31 p = 0.14	0.38 p > 0.5
Height Variation Contour [mm]	0.40 p > 0.5	0.55 p > 0.5
Mean RNFL Thickness [mm]	0.21 p > 0.5	0.22 p > 0.5

✔ Within normal limits  
⚠ Borderline  
✘ Outside normal limits

---

Pattern Deviation

GHT: Outside normal limits

Test:	24-2	24-2
Strategy:	ASTA Standard	ASTA Standard
FP Count:	3	3
FN Count:	3/3 (100.0%)	3/7 (42.9%)
Fixation Losses:	3 (0.7%)	0 (0.0%)

MN: \_\_\_\_\_

-7.24 dB p < 0.05%  
+5.66 dB p < 0.05%

● p > 5%  
● p < 5%

---

RNFL & Visual Field Combined OU Report

SPECTRALIS® & HEP

Patient: \_\_\_\_\_ DOB: \_\_\_\_\_ Sex: \_\_\_\_\_

Patient ID: \_\_\_\_\_ Comment: \_\_\_\_\_

Diagnosis: \_\_\_\_\_

**HEIDELBERG  
ENGINEERING**

**OD - Visual Field - 18.08.2010**

Stimulus: FDF Pattern: 24-2 Strategy: ASTA Standard

GHT Classification: Outside Normal Limits

Time span bet. RNFL and VF exams: 1 days

**OS - Visual Field - 18.08.2010**

Stimulus: FDF Pattern: 24-2 Strategy: ASTA Standard

GHT Classification: Outside Normal Limits

Time span bet. RNFL and VF exams: 1 days

**OD - RNFL Thickness - 18.08.2010**

IR 30°ART [HR]

Classification: Within Normal Limits

**OS - RNFL Thickness - 18.08.2010**

IR 30°ART [HR]

Classification: Within Normal Limits

Notes: \_\_\_\_\_

Date: 12.07.2011 Signature: \_\_\_\_\_

Software Version: HRT: 3.2/6133 HEP: 2.1.0.63


HRT & HEP –  
Unique Structure-Function Map

Image Courtesy Dr. Fritz Dannheim, Seevetal-Hittfeld, Germany and  
Dr. Victor Spears, Toronto, Canada.

SPECTRALIS & HEP –  
Combined RNFL & VF Report

# Precision

The SPECTRALIS enables the precise measurement of the peri-papillary retinal nerve fiber layer thickness. Comparing the measurement to a normative database allows a fast and reliable detection of local nerve fiber defects. The examination results are comprehensively documented and summarized.

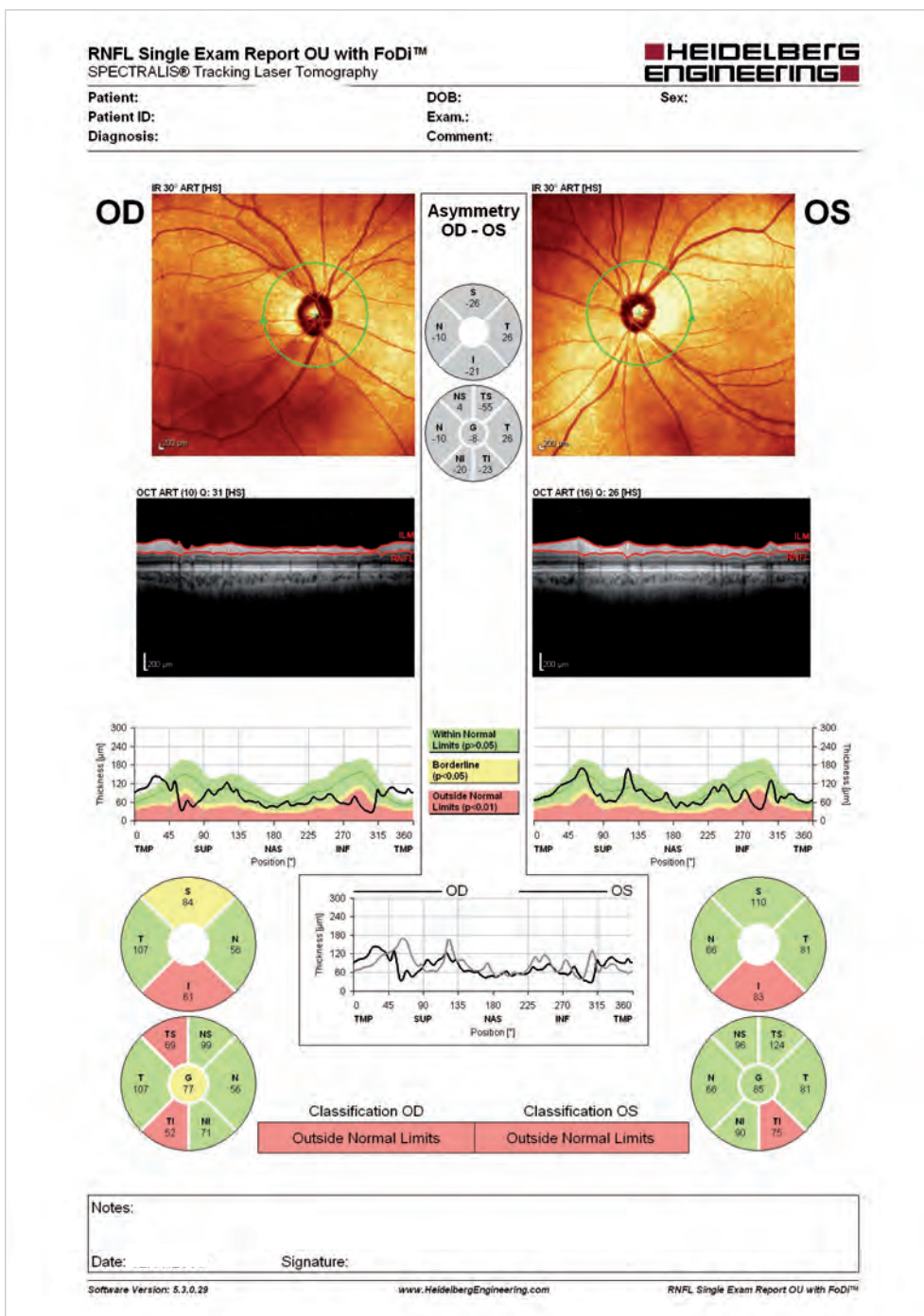


- ✚ Precise measurement of retinal nerve fiber layer thickness
- ✚ Evaluation of more than 1,500 data points on a single circle scan around the optic nerve head
- ✚ Automatic comparison to an age-adjusted normative database to classify a patient as being within or outside normal limits
- ✚ High sensitivity and specificity with TruTrack™ Active Eye Tracking and Fovea-to-Disc (FoDi™) alignment
- ✚ High-quality fundus images even under challenging conditions (e.g. cataract)



## Reproducibility

The reliable detection of change is critical and defines glaucoma in suspect patients. With the SPECTRALIS changes in nerve fiber layer and retinal thickness as small as  $1 \mu\text{m}$  can be identified.<sup>7,8</sup> Structured reports enable a comprehensive documentation of change over time.

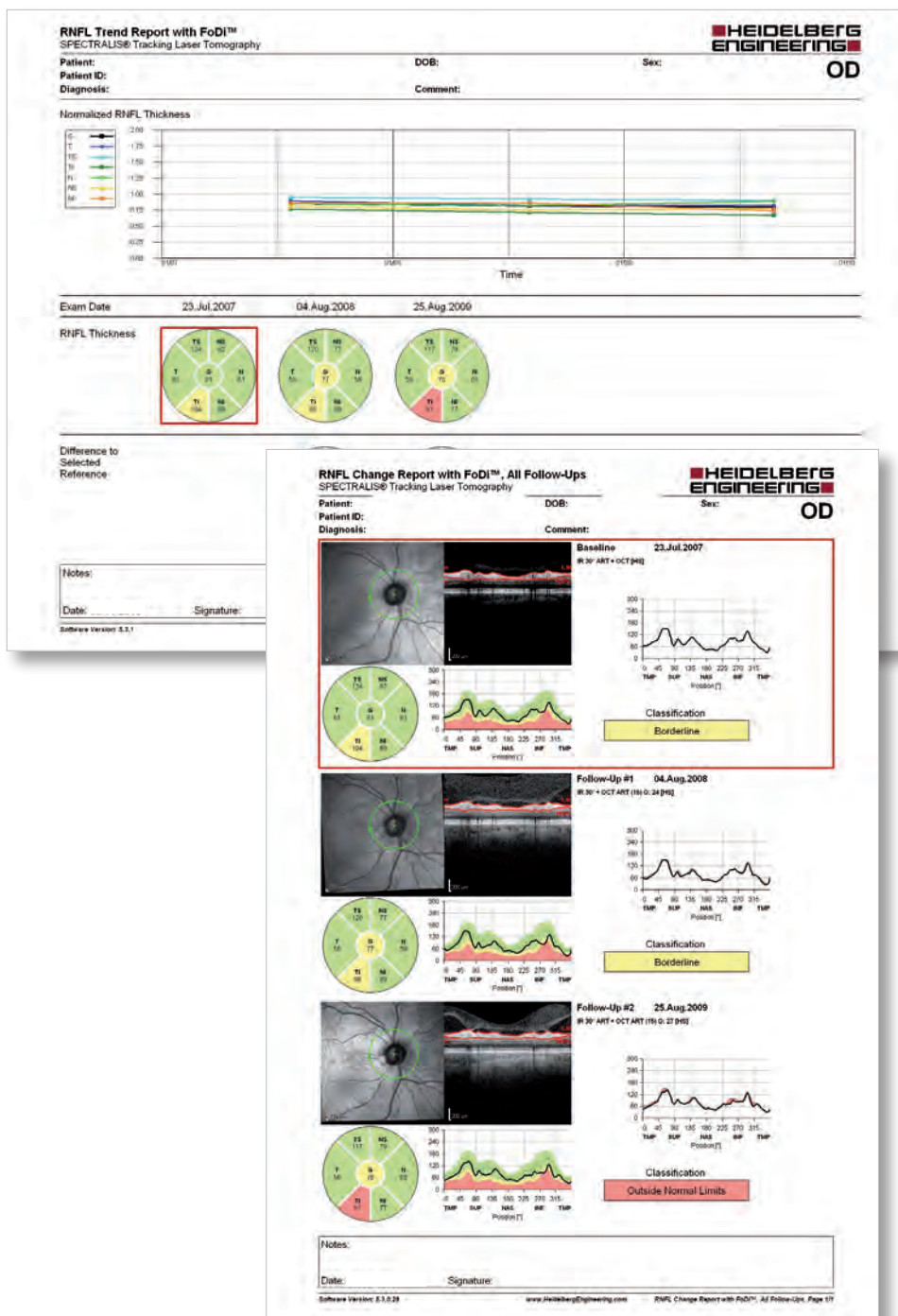

Reliable identification of changes as small as  $1 \mu\text{m}$

AutoRescan™ function for follow-up exams

Color-coding of local changes

Visualization of disease progression for easy interpretation

Comprehensive documentation of change over time


Fast localization of nerve fiber layer changes

<sup>7</sup> Wolf-Schnurrbusch et al. Macular thickness measurements in healthy eyes using six different optical coherence tomography instruments. Invest Ophthalmol Vis Sci. 2009; 50:3432-3437.


<sup>8</sup> Wu et al. Reproducibility of retinal nerve fiber layer thickness measurements using spectral domain optical coherence tomography. J Glaucoma. 2010 (Published online ahead of print).

## Innovative Concepts

The Posterior Pole Asymmetry Analysis combines mapping of retinal thickness over the entire posterior pole with an innovative asymmetry analysis. The full geographic extent of nerve fiber and ganglion cell layer defects can be visualized. The combination of nerve fiber layer thickness and retinal thickness provides an even more robust parameter for glaucoma detection and monitoring of change over time.

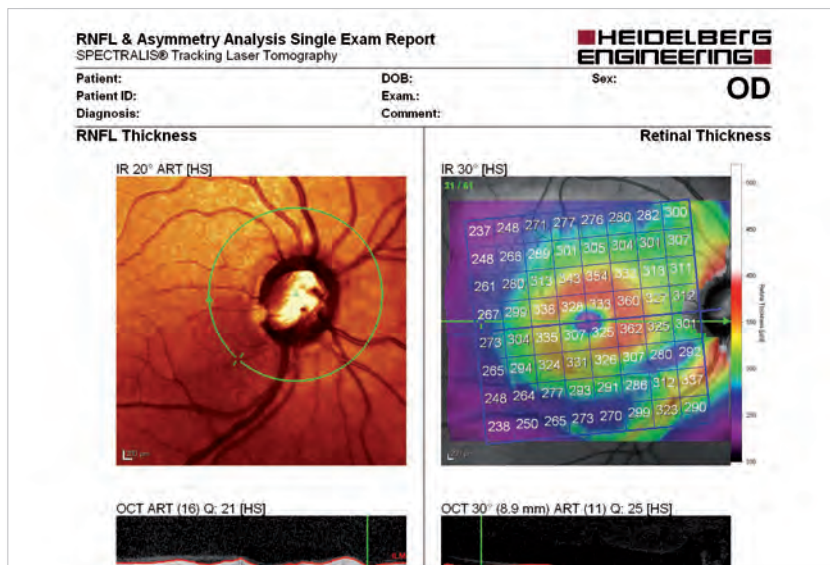


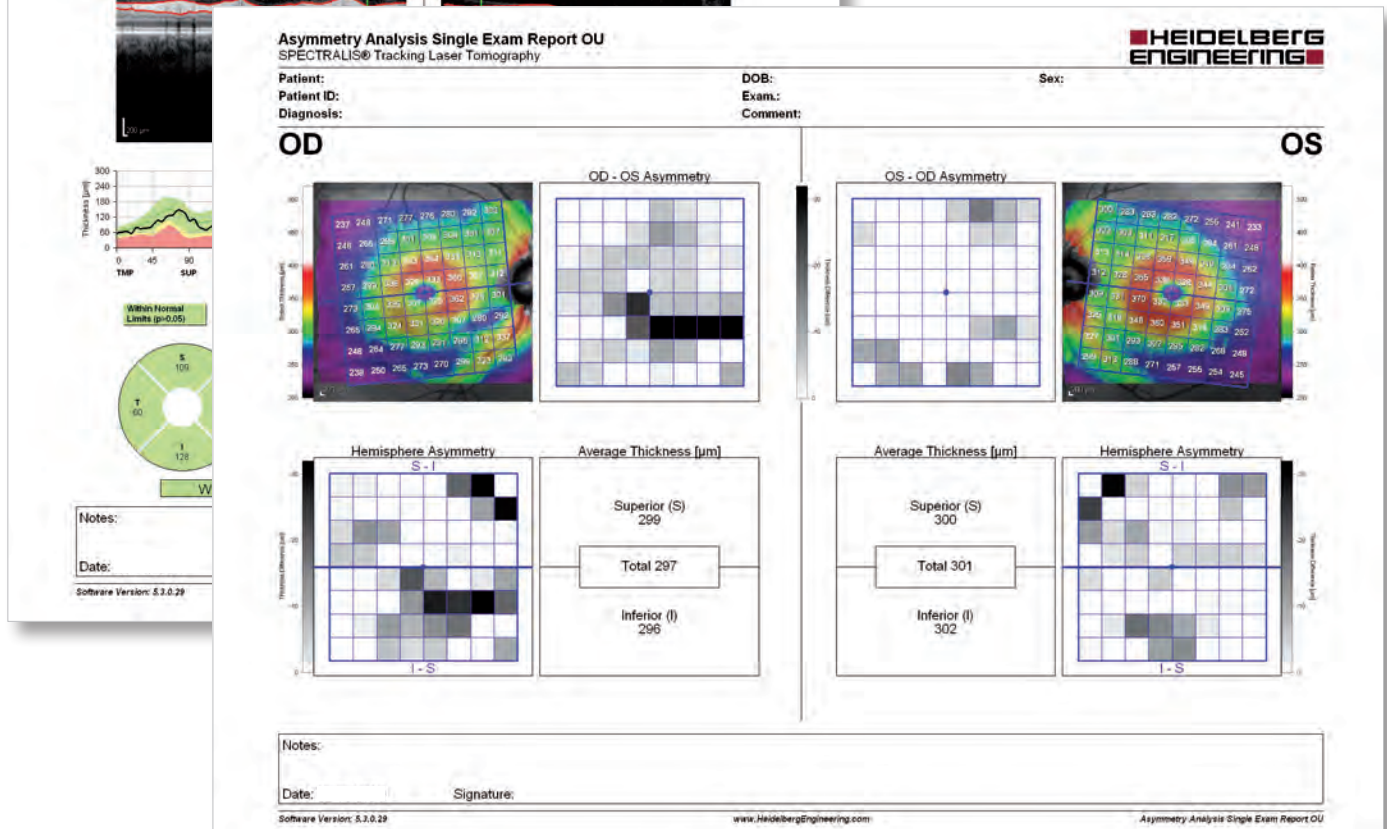


- ⚙ Evaluation of nerve fiber layer and retinal thickness in a single visit
- ⚙ Targeted combination of multiple parameters
- ⚙ High-resolution thickness maps with finely graded color scale
- ⚙ Asymmetry Analysis for early detection of glaucomatous abnormalities
- ⚙ Detailed reports



### Excellent Image Quality

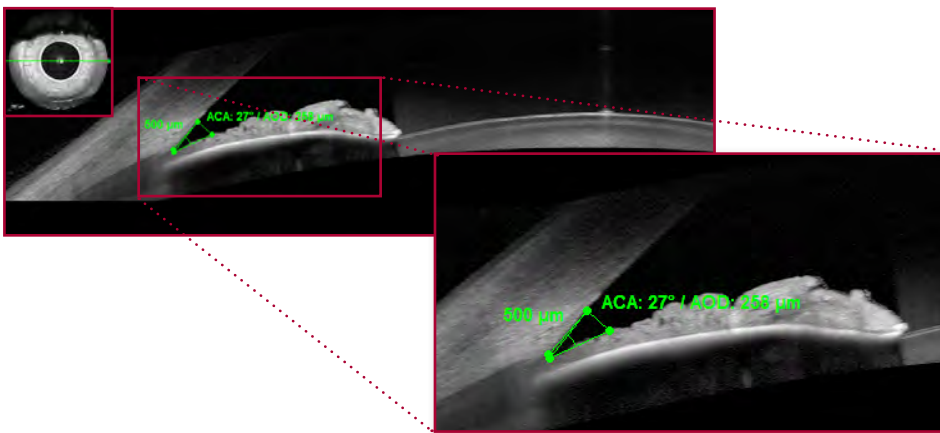Anterior Chamber Angle assessment is typically used in a comprehensive glaucoma diagnosis to exclude Angle Closure Glaucoma. The SPECTRALIS Anterior Segment Module offers image acquisition of the Cornea and the Sclera. Moreover, both anterior chamber angles can be visualized in a single SD-OCT scan for an efficient glaucoma diagnosis.



-  16 mm SD-OCT scan
-  High-resolution images acquired with TruTrack™ Active Eye Tracking and Heidelberg Noise Reduction
-  Manual measurement tools for Anterior Chamber Angle assessment

**Overview Report**  
SPECTRALIS® Tracking Laser Tomography

**HEIDELBERG ENGINEERING**

Patient:  
Patient ID:  
Diagnosis:

DOB:  
Exam.:  
Comment:

Sex:

**OD**

IR 30° ART + OCT 30° (16.7 mm) ART (60) Q: 45 EDI [HR]



IR 30° ART + OCT 30° (16.7 mm) ART (54) Q: 46 EDI [HR]



Notes:

Date: \_\_\_\_\_ Signature: \_\_\_\_\_



**Corporate Headquarters**

Heidelberg Engineering GmbH · Tiergartenstr. 15 · 69121 Heidelberg · Germany  
Tel. +49 6221 6463-0 · Fax +49 6221 646362

**UK Office**

Heidelberg Engineering Ltd. · Breakspear Park Suite F · Breakspear Way · Hemel Hempstead · Hertfordshire HP2 4TZ  
Tel. +44 1442 345 370 · Fax: +44 1442 345 001

**US Office**

Heidelberg Engineering, Inc. · 1808 Aston Avenue, Suite 130 · Carlsbad, CA 92008  
Tel. +1 760-536-7000 · Fax +1 760-536-7100

[www.HeidelbergEngineering.com](http://www.HeidelbergEngineering.com)